Instructions For Use FSG-Imprint-Foam

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1. Preparation

note the name of the customer and date the impression was taken on the foam box.
Open the box and with the Foam sides facing up at one. Lay a solid, level surface.

2. Heel

The patient's foot is centered on the foam sheet place and carefully on the Deposit foam surface. Well the check the perfect position of the foot. If the foot is not on the foam everywhere, correct now. Start with the By guiding your hand deep into the heel Press foam.

3. Metatarsus

Push the base M5 (midfoot) with the edge of the hand.


4. Big toe joint

Now press the big toe joint as well with the back of your hand in the foam. Relocate the power of your hand afterwards on your edge of your hand and use it to press the individual toes also.

5. Completion / repetition

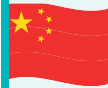Now take the foot vertically upwards from the foam. Repeat the the o.g. Processes with the other foot. Usability check the imprint, if necessary with a new one Repeat the pair of impression foams. Foam residues on the foot are made hygienic Reasons to remove immediately. FSG impression foam store separately with imprint.

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Formschaum Gronau GmbH